

User Manual

Q-LAS



A.R.C. Laser GmbH


CE 0123

Bessemerstr. 14
90411 Nürnberg

Tel.: +49 911 21779-0
Fax: +49 911 21779-99

info@arclaser.de
www.arclaser.de

Please contact for service:

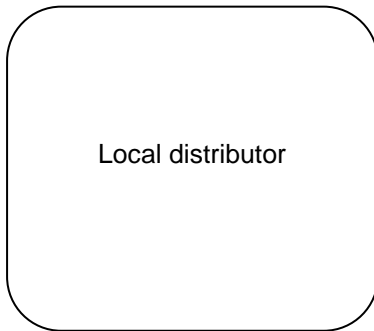

A.R.C. Laser GmbH Headquarters:

A.R.C. Laser GmbH

Bessemerstr. 14

D-90411 Nürnberg

Tel.: +49 911 21779-0

Fax: +49 911 21779-99

service@arclaser.de

Copyright © A.R.C. Laser GmbH. All rights reserved.

This document contains confidential and copyrighted information and may only be used with the devices described herein. Reproduction is prohibited without written permission of A.R.C. Laser GmbH.

1	Introduction.....	5
1.1	Marking and Symbols.....	5
1.2	Intended Purpose.....	6
1.2.1	General Purpose.....	6
1.2.2	Medical Purpose.....	6
1.3	General Contraindications.....	7
1.4	Theory and Technical Set-up.....	7
2	Transport and Storage	8
2.1	Shipping and Unpacking the Device	8
2.2	Return Shipment	9
3	Set-up and Installation	10
3.1	Installation Site	10
3.2	Room Requirement.....	10
3.2.1	Marking Access Points	10
3.2.2	Window Shielding	10
3.2.3	Reflecting Surfaces.....	11
3.3	Electrical Connection	11
4	Safety Information	12
4.1	General.....	12
4.2	Eyewear protection	13
4.3	Electrical Protection	13
4.4	Explosion and Fire Hazard.....	13
4.5	Protection against Undesired Radiation.....	13
4.6	NOHD Safety Distance	14
4.7	CE-Regulations	14
4.8	RoHS3-Regulations	14
4.9	Protecting Housing.....	14
4.10	Connectors and Switches	15
4.11	External Interlock Switch.....	16
4.12	Safety Shutter/Aiming Beam	16
4.13	Manual Reset	16
4.14	Reset by Power Failure.....	16
4.15	Labels and Markings	17
4.15.1	Type plate	17
4.15.2	Modification Label.....	17
4.15.3	Warning label laser.....	18
4.15.4	Serial number	18
4.15.5	Laser STOP	18
4.15.6	Follow instruction for use.....	18
4.16	Operating Conditions	19
4.17	Electromagnetic Compatibility	19
5	Advice for Users	20
5.1	Technical Instruction	20
5.2	Laser-Safety-Training.....	20
5.3	Medical Instruction	20
5.4	Device Parts and Accessories	21
5.4.1	Basic equipment	21
5.4.2	Optional accessories	21
5.4.3	Instrument Table.....	22
5.4.4	Foot Switch (optional)	23
5.4.5	Lowering prism (optional)	23
6	Operation	24
6.1	Starting the Laser	24
6.2	Display.....	25
6.2.1	Energy Selection.....	25
6.2.2	Burst Mode (Number of pulses = pulses per shot)	26

6.2.3	Information LED Energy Level.....	26
6.2.4	Reset Counter.....	26
6.2.5	Laser STANDBY / READY	26
6.2.6	Counter and cumulated applied energy.....	27
6.2.7	DuoSpot Aiming Beam	27
6.2.8	Head-Up Display.....	27
6.3	Adjustments at the laser head	28
6.4	Trigger Button at the Joystick.....	28
6.4.1	Standard version with rotating prism	28
6.4.2	Version with lowering prism	28
6.5	Treatment.....	29
6.6	Switch Off.....	29
6.7	Laser-Stop.....	29
7	Specifications	30
7.1	General.....	30
7.2	Laser Data.....	30
7.3	Electrical Connection Data.....	30
7.4	Classification	30
8	Service.....	31
8.1	Introduction.....	31
8.2	Safety Check (STK)	31
8.3	Care and Maintenance by the User	32
8.3.1	Inspection and cleaning of external optics.....	32
8.3.2	Inspection and cleaning of internal Slit lamp optics.....	33
8.4	Slit lamp.....	33
8.4.1	Dust Cover	33
8.4.2	Test the Slit lamp focus	34
8.5	Troubleshooting	35
8.6	Disposal.....	36
9	Customer Service	37
9.1	Warranty Information.....	37
9.2	Warranty, Shipment, Packing	37
9.3	Sales and Service Information	37
10	Guidelines and Manufacturers Declaration	38
10.1	Electromagnetic Emissions	38
10.2	Electromagnetic Immunity (1)	39
10.3	Electromagnetic Immunity (2)	41
10.4	Recommended separation distances between portable and mobile RF telecommunications equipment and the laser	42

1 Introduction

We thank you for having chosen our Q-LAS and wish you every success in using this versatile laser system.

The Q-LAS laser generates a laser beam of high intensity, which can cause injuries if handled improperly. The user manual should therefore be read carefully before starting up the device. If you have any further questions regarding the safety or the use of the device or regarding the laser and laser radiation, please contact A.R.C. Laser GmbH.

1.1 Marking and Symbols



The "General warning sign" (ISO 7010 W001) indicates general risk. It encourages the user to take care regarding the hazard specified by the supplementary sign. Before carrying out any further work on such marked parts, the operating instructions should be consulted or the service of A.R.C. Laser GmbH should be contacted.



The sign "Warning; Laser beam" (ISO 7010 W004) indicates laser hazards. It encourages the user to take care to avoid exposure to a laser beam. The laser beam of this device is not visible to the naked eye nor through protecting goggles.

1.2 Intended Purpose

1.2.1 General Purpose

Tissue disruption by plasma shockwave in the human eye.

1.2.2 Medical Purpose

The Q-LAS is a medical laser device intended to be used in ophthalmic applications particularly for the treatment of glaucoma (closed angle glaucoma) and posterior cataract opacification (PCO) patients. The Q-switched Nd:YAG laser delivers very short but high energetic laser pulses. The pulse length is in the nanosecond range with adjustable pulse energies. The short, high-energy pulses generate a plasma in the focal plane, which provides an optical breakthrough. Regarding to the medical indication (closed angle glaucoma or PCO) the Q-Las provides the treatments of iridotomy or capsulotomy.

Medical Application	Indication	Side and adverse effects	Contraindications
Angle closure glaucoma (Iridotomy)	Prevention or treatment of high intraocular pressure related to angle closure glaucoma	<ul style="list-style-type: none"> • IOP rise (non-permanent) • inflammation • bleeding • change of endothelial cell density • cystoid macular edema • retinal complications • scotoma • discomfort (mild/ moderate) • pain (mild) 	<ul style="list-style-type: none"> • open angle glaucoma • congenital glaucoma • inflammatory/ uveitic glaucoma • acute intraocular inflammation • acute reduced optical transparency of the cornea
Posterior cataract opacification (PCO) (Capsulotomy)	reduced visual acuity related to PCO	<ul style="list-style-type: none"> • inflammation • change of endothelial cell density • cystoid macular edema • retinal detachment • IOP rise • bleeding • lens damage/ pitting lens 	<ul style="list-style-type: none"> • acute intraocular inflammation • acute reduced optical transparency of the cornea

ATTENTION

The device may only be operated by specially trained personnel who are experts in the medical effects and possible dangers of the device. You must have the necessary skills to use the laser in accordance with this instruction manual.

When not in use, the device should be protected against unauthorized use.

1.3 General Contraindications

The device may only be used for iridotomy and capsulotomy. Any other treatments are not intended with Q-LAS.

1.4 Theory and Technical Set-up

The Q-LAS is a Nd:YAG laser with a wavelength of 1064 nm. The radiation in the laser head is generated with the help of a diode. The laser beam that leaves the device has a small diameter; this is less than 10 µm at the place of its effect in the patient's eye.

The wavelength of the emitted laser light depends on the laser medium. In case of the Q-LAS, infrared laser radiation is generated. Attenuating optics allow to adjust the output energy. Internal control mechanisms compare the actual and the target state for the selected amount of energy.


With the Q-LAS, the laser beam path is overlaid with the observation beam path of the slit lamp. With the help of a target beam, the desired treatment position can be precisely determined in advance. Furthermore, the Q-LAS has a retractable prism.

The focus of the laser can be shifted in 2 steps using a dial on the laser unit.

2 Transport and Storage

We at A.R.C. Laser GmbH will make sure that the device is packed and transported with the greatest possible care.

Before unpacking the laser, please check the packaging for damage and report any damage immediately to the shipping agent and A.R.C. Laser GmbH. Only remove the packaging in the presence of a representative of the carrier. Make a list of the damaged parts and have this list signed by the courier.

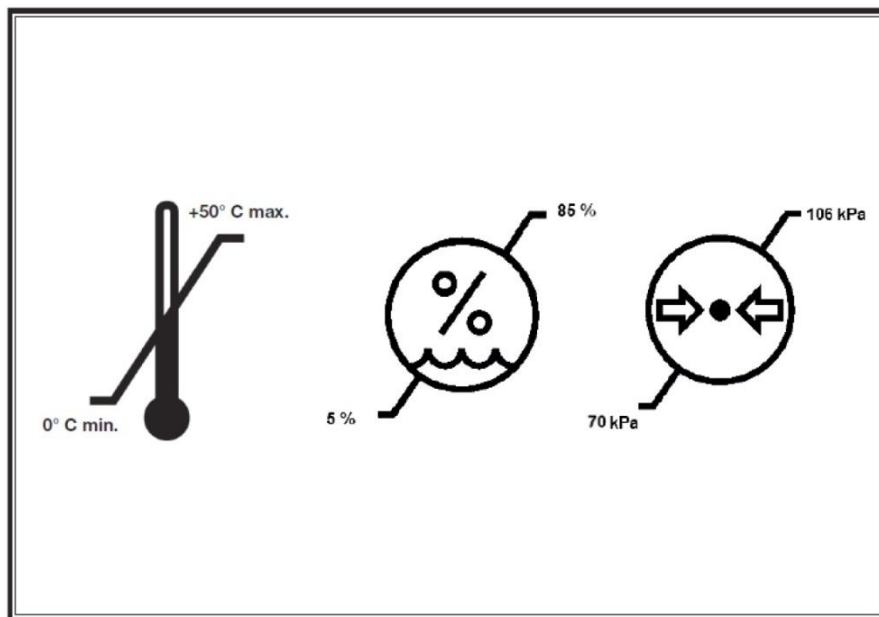

The device must be transported at temperatures between -10°C and +50°C. The air pressure during transport must be between 1080 hPa and 750 hPa. During storage, a temperature range of -10°C to 50°C must be maintained. The environment/air must be dry and clean. The relative humidity during transport and storage must be between 10% and 90%.

ATTENTION

If the laser is transported or stored at a temperature below 5 °C, it can be damaged when starting. Unpack the laser and leave it for at least 12 hours at normal room temperature so that the system reaches room temperature.


2.1 Shipping and Unpacking the Device

Unpacking and installation of the device must be carried out by an authorized technician or authorized representative of A.R.C. Laser GmbH. After the installation and the correct connection of the device, the technician will put the device into operation and explain the function of the device. All functions and safety procedures are discussed and demonstrated here.

2.2 Return Shipment

The device will be shipped in a specially developed transport packaging. Please keep this. If you ship the device back to A.R.C. Laser GmbH, be sure to use the intended transport packaging.

The transport packaging is designed so that the device fits snugly into the foam parts.

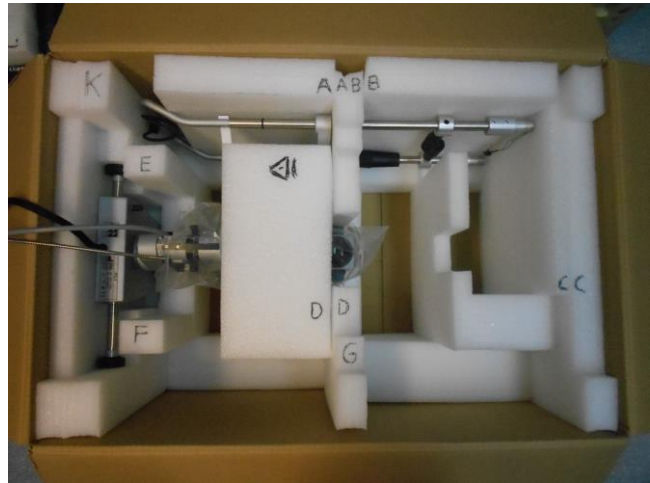

All accessories are stowed in an extra padding inside the transport-packaging.

If all foam parts are in the right place and the box is closed, the device is ready for shipping.

3 Set-up and Installation

3.1 Installation Site

Ensure to provide a suitable location for the device before it is set up.

The Q-LAS should be operated in an easily accessible place. The laser should not be operated near a heater, since the device only works when the ambient temperature is not higher than 45°C. Higher temperatures can cause the device to switch off due to overheating. A place of installation in direct sunlight can cause an early shutdown and must be avoided in any case.

If the temperature is too low (below 10°C), the device cannot be started to prevent condensation on the internal optics; this could result in permanent damage to the laser.

- The device should be set up so that the laser beam is not aimed directly at a door, window or reflective material.
- The device should be operated in a dust-free room. There should be no carpets on the floor or the walls.
- When the device is not in use, it should be covered to prevent contamination. For this purpose, a suitable dust cover is included.

The wiring must be installed in a way that no tripping hazards or other hazards occur.

All control elements must also be free and easily accessible. The air humidity is monitored internally and must be below 75%.

3.2 Room Requirement

The legislation imposes the following requirements to any room in which a class 3B laser according to IEC 60825-1 is operated.

3.2.1 Marking Access Points

All entrances must be clearly marked so that unintentional or unauthorized entry, which can lead to a hazard, is avoided.

- Laser warning signs (triangle with laser symbol) and the wavelength labeling must be attached to each access door.
- A warning lamp must be attached above each access door. This must always light up when the laser is in operation. Unintentional entry into the room without safety glasses is avoided.
- The laser safety goggles must be easily accessible at the entrance.

3.2.2 Window Shielding

It must be ensured that no laser radiation can escape the room. In particular, the windows must be secured with suitable materials. If you have any questions, your contact person at A.R.C. Laser GmbH is available to assist you.

3.2.3 Reflecting Surfaces

To avoid danger from reflected direct and scattered radiation, there must be no reflective surfaces in the room.

These can be:

- Mirrors
- Pictures behind glass
- Chrome surfaces
- Windows

These surfaces must be removed or suspended or matted. In the area of laser use, only matted, non-reflective and non-flammable instruments and materials should be used.

3.3 Electrical Connection

The laser can be operated from 100 V - 240 V (50 Hz/60 Hz) AC voltage.

The power connection is made via the table column and the supplied IEC connector.

Make sure that the plug is accessible at any time so that the laser can be disconnected from the mains after use.

4 Safety Information

4.1 General

The Q-LAS laser is a precision instrument for medical applications. The system has been carefully developed and tested through intensive testing before shipping. In order to offer you and the operating personnel all possible protection, we recommend that you read this section of the operating instructions carefully.

According to EN 60601-2-22, respectively EN 60825-1 the Q-LAS is categorized into laser class 3B.

Class 3B lasers are high-energy lasers. To ensure safe and trouble-free working with the device, special arrangements must be made. It is particularly important to protect the eyes and skin of the operator, the patient and third parties. Laser safety goggles must be used for eye protection.

The following explanations are not exhaustive. All users of laser devices should enclose applicable legal regulations and provisions with the device and inform the staff accordingly.

If the device is operated outside Germany, provisions of the American National Standard Office ANSI Z136.3-2018 "American National Standard for the Safe Use of Lasers in Health Care Facilities" and ANSI Z136.1-2014 "American National Standard for the Use of Lasers" must be respected.

This manual is limited to the operation, maintenance and control of the device. The manual is not a guide for the treatment of diseases that can be treated by laser.

Regarding the additionally supplied devices, such as slit lamps or instrument tables, the safety, operating and maintenance instructions in the corresponding manuals must be followed. Regarding the instrument table, reference is only made here to the danger from improper use of the table as a seat or storage area. When operating the height adjustment, make sure that no one can be harmed.

A.R.C. Laser GmbH cannot be held responsible for damage or damage resulting from improper use.

The warranty of the device expires if the laser has been opened (even partially), modified or repaired by unqualified personnel.

4.2 Eyewear protection

As a safety measure against direct or indirect laser radiation, it is necessary that everyone in the room wear laser safety goggles. The operating doctor is protected from radiation by the laser protection integrated (eye safety filter) into the slit lamp. Appropriate eye protection of the patient must be guaranteed.

When using the laser, only those laser safety goggles that are designed for the wavelength of 1064 nm may be used. The safety goggles should bear a CE mark as well as the appropriate protection class.

The following laser safety goggles from A.R.C. Laser GmbH are suitable for the Q-LAS:

- AS01033: Protection class D LB 6 + IR LB8 + M LB8Y / Optical density OD 8+

Safety goggles can be worn even when the person wears glasses. Scattered radiation that does not reach the eye directly from the front may present risks resulting from internal reflections on the glasses. Therefore, we recommend wearing safety goggles that also ensure full side protection. The safety glasses offered by A.R.C. are suitable as over-glasses for spectacle wearers and also provide lateral protection.

ATTENTION

Never look directly at the laser beam or the light reflected by the laser beam. This may cause serious damage to your eye.

4.3 Electrical Protection

Never remove the housing parts of the laser. Only authorized A.R.C. Laser GmbH service personnel is allowed to modify the device. If the mains have to be replaced, the replacement must be carried out by authorized service personnel.

The room in which the laser is operated should be kept dry. After cleaning the room, make sure the floor is dry before starting up the device.

ATTENTION

Never use the laser if you notice any visible defect on the device.
Never use the device if wires or supply lines are exposed.

4.4 Explosion and Fire Hazard

Never work with the laser near flammable anesthetics, flammable solutions or other flammable materials. Flammable plastic parts or paper parts should also be removed from the working area of the laser in particular. There is a fire or explosion hazard if the laser is used in the presence of flammable materials, solutions, or gases or an oxygen-enriched environment.

4.5 Protection against Undesired Radiation

If a footswitch is used, it should never be outside the area of the attending doctor. It is forbidden for anyone other than the attending doctor to trigger the footswitch.

Especially in operating theaters where there are multiple footswitches, it is important to ensure that the laser footswitch is close to the laser device. Also, the hand trigger of the slit lamp may only be used by the attending surgeon.

4.6 NOHD Safety Distance

The NOHD (Nominal Ocular Hazard Distance) is the distance at which the irradiance is equal to the exposure limit value of the cornea of the eye. The NOHD identifies the danger area within which there is a risk of damage to the health of the eye if you look directly into the laser beam without protection.

The NOHD is calculated according to EN 60825-1 and taking into account the permitted power fluctuations (+/- 20%) according to EN 60601-2-22 using the following equation:

$$NOHD = \frac{\sqrt{\frac{4P}{MZB * \pi}} - Diameter\ Beam\ Bundle}{Beam\ Divergence}$$

For the Q-LAS the following NOHD applies:

Wavelength λ:	1064 nm
NOHD:	3.6 m
Beam divergence $\alpha_{(full\ angle)}$:	16°
Permitted Maximum Radiation:	2 x 10⁻² Jm⁻²

4.7 CE-Regulations

The laser system Q-LAS has been approved by the notified body in accordance with the European Regulation 2017/745/EU for medical devices. Accordingly, the device bears the CE mark **CE 0123**.

The device has been checked for electrical and mechanical safety. All parts we use comply with the CE regulations or have been tested for approval or suitability by the notified body. Additional devices that you attach to the device require approval from an official test center. Changes to the device or interventions on your part will void the approval and warranty. The test approval number is included with the device.

4.8 RoHS3-Regulations

Our company operates worldwide and regards the protection of the environment and natural resources as an entrepreneurial obligation. Based on individual tests, A.R.C. Laser GmbH confirms that, to the best of our knowledge, our products do not contain any substances in concentrations whose placing on the market is prohibited according to the applicable requirements of Directive 2015/863/EU (RoHS3).

4.9 Protecting Housing

The Q-LAS laser system has a protective housing. It prevents radiation of the laser from escaping and protects users from touching energized parts. This housing must not be removed.

4.10 Connectors and Switches



Main switch

Thermal fuses



- Door interlock
- Warning light
- Foot switch/
hand release
- Connection laser device
- Power supply
- Connection of the laser
shot LED



Laser stop and key switch

The main switch under the table turns the entire system on and off. The key switch is used to turn the laser on. The stop button is used to switch off the laser in emergencies.


To prevent a defect in the unit - caused by the power supply unit - the Q-LAS has two thermal fuses under the table. When the thermal fuses are tripped, they are pushed down and the red LED on the main switch fades out. As a result, no more voltage is passed on to the Q-LAS and you cannot continue treatment. In this case, please contact the service department of A.R.C. Laser GmbH or their local representative.


The sockets of the device are located on the bottom of the instrument table.

The red marked connector is the interlock connector ○. You can replace it with a door interlock switch (see chapter 4.11).

You can connect a warning light to the orange marked connector ○. This indicates to third parties whether the laser is currently in operation. The switch and wiring must be designed for at least 12 V and 100 mA.

The manual release or the foot switch is connected to the blue marked connection .

The Q-LAS is connected to the table at the green marked connection .

The yellow marked connection is the power supply for the Q-LAS .

The purple connector  is used for the laser shot LED.

During the installation of the device, the authorized technician or representative of A.R.C. Laser GmbH will connect all needed wires. The user does not need to change any connections. If problems arise, please contact the A.R.C. Laser GmbH or the local representative.

4.11 External Interlock Switch

A door interlock switch is required by the accident prevention regulations. The device is equipped with an interlock connector as standard, which can be replaced by a door interlock switch. The laser switches off when the door is opened. In addition, the laser cannot be switched to the READY mode when the door is open. An error message appears on the control panel. When the door is closed, the error message disappears and the laser can be switched READY again.

When installing a door interlock switch, note the following:

The switch and wiring must be provided for at least 12 V and 20 mA. The wires should end with a standard male connector. The choice of polarity is expedient in both variants. The socket for connecting the door switch can be found underneath the table. Insert the door switch connector there.

Make sure the socket sticks firmly to prevent unexpected system interlock problems.

4.12 Safety Shutter/Aiming Beam

The Q-LAS has an internal safety shutter. This safety switch enables the release of the laser beam mechanically. As soon as the laser has passed its internal tests and calibrations, the laser can be set to READY mode by pressing the READY button. As soon as the laser is set to READY mode, the aiming beam will appear. The aiming beam is a low-level laser with very low power.

4.13 Manual Reset

If an error occurs, the system switches into STANDBY mode. You should then switch the laser off and on again. The restart should fix the error due to the automatic recalibration. Please refer to chapter 8.5 or a description of possible errors.

If the error still occurs, this can only be remedied by trained personnel. Please contact the service of A.R.C. Laser GmbH in this case.


4.14 Reset by Power Failure

If the device is disconnected from the power grid accidentally - e.g. due to a power failure - it usually restarts automatically. The system recalibrates and deletes all unsaved settings. If there is still an error, it can only be corrected by trained personnel. Please contact A.R.C. Laser GmbH service department should this occur.

4.15 Labels and Markings

The Q-LAS is provided with various warnings in accordance with European and worldwide guidelines. This is to prevent the laser user from being exposed to laser radiation due to carelessness. Following labels are attached on the Q-LAS:

4.15.1 Type plate

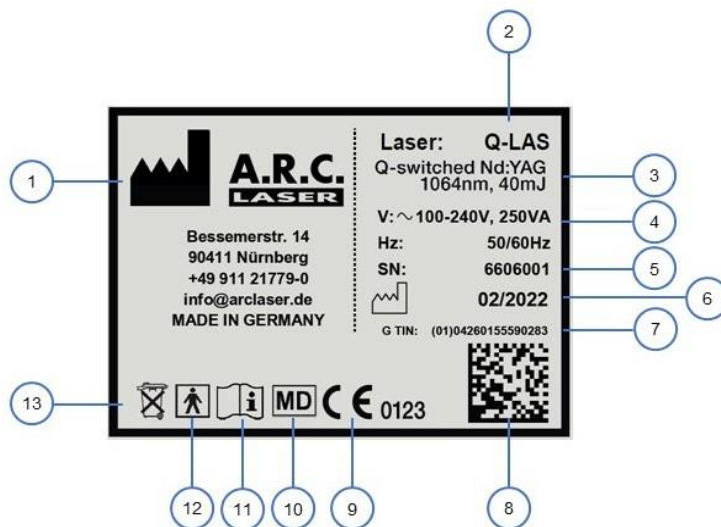

- 1) Manufacturer
- 2) Laser name
- 3) Laser specifications
- 4) Electrical connection data
- 5) Serial number (SN)
- 6) Production date
- 7) GTIN (Global Trade Item Number)
- 8) UDI (Unique Device Identification = GTIN + SN + production date)
- 9) CE Mark
- 10) Medical Device
- 11) Follow instructions for use
- 12) Application part type BF
- 13) Do not dispose in household trash

4.15.2 Modification Label

MODIFIKATION					
1	2	3	4	5	6
7	8	9	10	11	12
13	14	15	16	17	18
19	20	21	22	23	24

Shows the current device status based on the marked modification numbers.

4.15.3 Warning label laser

Laser Aperture



The “Laser Aperture” label marks the laser output. The other labels provide information about the type and intensity of the laser system.



4.15.4 Serial number

A.R.C. LASER PCL5 SN: 12010593	A.R.C. LASER Laserbox Q-LAS SN: 68001
Part of Q-LAS, SN: 6606001	Part of Q-LAS, SN: 6606001
A.R.C. LASER U-Table SN: 0058	A.R.C. LASER Q-LAS SN: 6606001
Part of Q-LAS, SN: 6606001	

The Q-LAS consists of several components, each with its own serial number.

4.15.5 Laser STOP



This symbol indicates where the laser STOP button is located. (see chapter 6.5)

4.15.6 Follow instruction for use



On the device, you will again be reminded that the instructions for use must be respected.

4.16 Operating Conditions

The medical laser Q-LAS is not suitable for use in connection with combustible gas mixtures of all kinds.

The device is not approved for operation at altitudes above 2,000 m above sea level. The air pressure must be between 1080 hPa and 750 hPa.

The following environmental conditions must be met:

- Ambient temperature: 15°C to 45°C
- Relative humidity: < 75%, non-condensing

4.17 Electromagnetic Compatibility

The Q-LAS laser system meets the EMC requirements according to EN 60601-1-2. Guidelines and the manufacturer's declaration are described in chapter 10.

ATTENTION

Avoid using this device next to other devices or with other devices in a stacked form as this could result in improper operation. If such use is still necessary, this device and the other devices should be monitored to ensure that they are operating properly.

5 Advice for Users

5.1 Technical Instruction

During the installation, the A.R.C. Laser employee or local authorized representative will give a briefing. It contains the technical use of the device and essential security-relevant points. During the briefing, every person working in the vicinity of the laser should be present.

5.2 Laser-Safety-Training

The Q-LAS laser is designed for medical users. It may only be used by personnel who have been instructed in its operation.

A.R.C. Laser GmbH recommends attending additional seminars to train working with different laser systems.

The A.R.C. Laser employee or authorized representative who installs the device will also offer additional training for accompanying personnel. Subject of this briefing are the usage of laser safety goggles and laser safety in particular.

A.R.C. Laser GmbH has a list of recommended courses as well as laser safety courses. These can be obtained from us at any time.

5.3 Medical Instruction

During the instruction of the device, the A.R.C. Laser employee or authorized representative will only give a general overview of the medial applications.

If necessary, there is the possibility to take part in a training course with an experienced physician. Please contact your local A.R.C. Laser authorized representative, or contact us directly.

5.4 Device Parts and Accessories

5.4.1 Basic equipment

Additionally to the device itself, the basic parts consists of functional parts that are necessary for the operation of the device, as well as additional equipment for safety and functionality.

Necessary parts for operation:

Part	Description	Article number
Mains cable*	Type E/F (EU)	KB13003
	Type G (UK)	KB13004
	Type B (US)	KB13005
Door interlock	Door interlock plug, angular	KB03005
Chin rest paper**	Chin rest paper for the patient	SL01115
User Manual	User manual Q-LAS	

*The type of the mains cable depends on the destination country. One cable is included in the standard scope of delivery.

**Only chin rest paper provided by A.R.C. laser is allowed to be used with the Cobra. When the enclosed chin rest paper is used up, refills can be obtained from A.R.C. laser.

Additional equipment included:

The packaging equipment consists of:

Part	Description	Article number
Transport carton	Shipping box incl. distance foam	VP06028
Transport box***	Wooden shipping box	VP01081

***Lasers that are not delivered via freight transportation may not include a wooden shipping box.

other equipment:

Part	Description	Article number
Dust protection cover	Dust protection cover for the device	ZU01046
Cover for toothed racks	Cover for toothed racks to adjust the slit lamp position	ME06544
Spare LED Modul	Spare module for the LED	BG06565
Focus rod	Focus rod for adjustment of binoculars	ME06883

5.4.2 Optional accessories

Optional accessories are marketed separately. The products listed can, but do not have to, be purchased from A.R.C. laser.

A.R.C. laser offers **laser safety goggles** for eye safety of additional people in the operating room.

Part	Description	Article-No.
Laser safety glasses	Protection class D LB 6 + IR LB8 + M LB8Y / Optical density OD 8+, 1064 nm	AS01033

Contact glasses* required for SLT, Iridectomy or capsulotomy procedures can also be obtained from A.R.C. laser:

Application	Description**	Article-No.
Iridectomy	Ocular Abraham Iridectomy YAG Laser Lens	WE01243
Capsulotomy	Ocular Abraham Capsulotomy YAG Laser Lens	WE01244

*All offered contact glass forms are well established for their respective treatment. For further information about the medical treatments or the recommended devices consult the application manual or your medical device consultant. When buying contact glasses for use with Q-LAS, make sure that they are suitable for laser treatment.

**The listed designations are the trade names by ocular instruments. They may differ for different providers.

For the personal comfort of the physician and the stabilization of the arm while holding a contact glass, any comfortable **elbow rest** can be used. A.R.C. offers a stackable elbow rest made of foam. (ZU11056, diameter 116mm, height 28mm).

ATTENTION

Only spare parts and applicators which are approved by A.R.C. Laser GmbH are to be used with the device. Accessories that have not been approved can significantly impair the safety and reliability of the device.

The use of accessories, transduce and wires other than those which the A.R.C. Laser GmbH has determined or provided, may result in increased electromagnetic interference or reduced electromagnetic immunity of the device and lead to incorrect operation.

5.4.3 Instrument Table

The instrument table on which the slit lamp is mounted is height adjustable. This allows you to adjust the optimal working height comfortably for each patient. Please note that the table can only take an additional load of 2.0 kg. This means that no objects heavier than 2.0 kg in total may be placed on the table. Also, no persons may sit on the table or support themselves with their entire body weight on it.

If the table is tilted, you can calibrate it. To calibrate the table, move it all the way down. Press and hold the button (down arrow) until it comes to a complete stop. After calibration, the table is horizontal again.

5.4.4 Foot Switch (optional)

An optional footswitch is available for your Q-LAS. This enables to emit laser radiation via the footswitch.

5.4.5 Lowering prism (optional)

A lowering prism can also be installed in the Q-LAS. This enables treatment by lowering the prism, without swinging out the lightning arm. A rotational prism is installed as standard.

6 Operation

This part of the manual only describes the technical application of the device, without going into the medical application.

Settings and adjustments should only be made in accordance with the operating instructions. Changes or settings that are not described in this manual can lead to malfunctions.

The laser is in the READY mode during treatment. If you interrupt, the laser must be returned to STANDBY mode. The device must always be switched off when unattended to prevent operation by an unauthorized third.

ATTENTION

Since the aiming beam takes the same way through the laser transmission system as the working beam, controlling its appearance is a good way to check the integrity of the laser transmission system. The aiming beam will appear as soon as the device is set to READY mode. If the aiming beam does not appear at the distal end of the laser transmission system, its intensity is weak or if it looks diffuse, the laser transmission system might be damaged.

ATTENTION

The use of the controls or adjustments for performing procedures other than those specified herein may result in hazardous laser radiation exposure


6.1 Starting the Laser

- Make sure the power cord is plugged in. The power cable must not be installed under tension.
- Make sure the door contact switch is connected.
- Check that the Laser-Stop button is not pressed.
If it is pressed, pull it up.
- Are third parties present who need safety goggles?
- Activate the main switch at the table.
- Turn the key switch.

ATTENTION

The Q-Las should not be operated for longer than 5 hours at a time. Restart the laser after 5 hours at the latest so that it can carry out the system check during the start routine. This ensures that undetected errors cannot occur.

6.2 Display

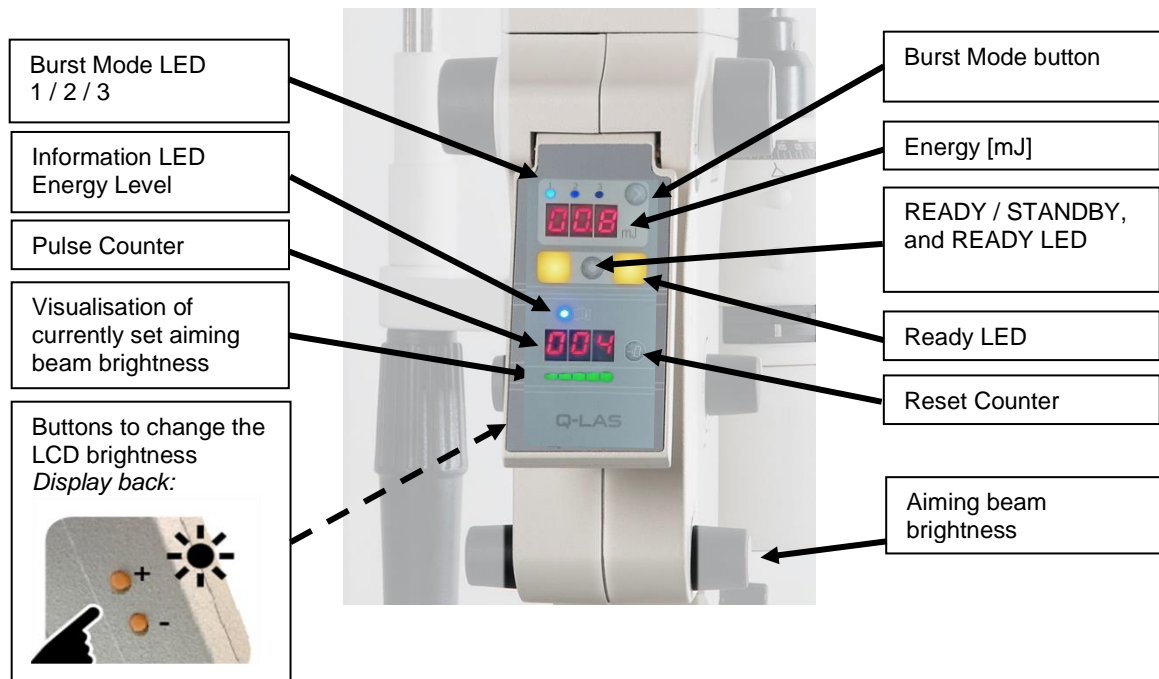

The display shows different information depending on the state of the unit.

When the Q-LAS is started by turning the key switch, the following sequence occurs:

- Warm-up phase: the system counts from 100 to 0. If the system has not yet reached the required operating temperature after the warm-up phase (approx. 60 seconds), it counts again from 100 to 0. A total of three warm-up phases are possible. After the warm-up phase, the device is at operating temperature (max. 180 s).
- The LEDs start to blink - - -
- After successful system check, the energy display shows the currently set value.
- If an error message is displayed, you can find out how to correct it in chapter 8.5.

ATTENTION

Please do not change any parameters during system start.


6.2.1 Energy Selection

Use the energy selection wheel to change the values. The available energy values can be continuously adjusted from 0.5 mJ to 10 mJ (BurstMode 1), 20mJ (BurstMode 2) up to a maximum of 30mJ (BurstMode 3).

6.2.2 Burst Mode (Number of pulses = pulses per shot)

The number of pulses can be varied between 1, 2 or 3. The total energy output (mJ) changes depending on the number of pulses.

By changing the burst mode, the total energy changes as follows:

Energy (mJ) according to burst mode (1, 2 or 3 pulses)		
1 Pulse	0.5 mJ to 10 mJ in 0.1 mJ steps	
2 Pulses	1.0 mJ to 20 mJ in 0.2 mJ steps	
3 Pulses	1.5 mJ to 30 mJ in 0.3 mJ steps	

6.2.3 Information LED Energy Level

This information LED is off as long as the laser energy is exposed within the standard tolerance range of +/- 20 %. As soon as the tolerance of the exposed laser energy is between + 20 % up to + 30 % the information LED lights up in blue and in addition an acoustic signal is emitted to inform the user that the laser is working in the extended tolerance range.



In case that the energy exceeds this extended tolerance range, the laser goes into the test mode. After the test mode has been successfully completed, the treatment can be continued. In case the test mode is not successfully completed a restart of the device is needed.

6.2.4 Reset Counter

To reset the counter to 0, press the Reset-button on the remote control: 000 is displayed again.



6.2.5 Laser STANDBY / READY

This button switches from STANDBY to READY and back. There is a safety break of 2 seconds programmed before the laser fires the first shot. During this safety break, the yellow LEDs flashes.

Yellow LED **lights up**: Laser READY

Yellow LED **does not light up**: Laser STANDBY



ATTENTION

If the laser is in READY mode and is not used for more than 2 minutes, the laser automatically switches to STANDBY mode.

6.2.6 Counter and cumulated applied energy

By default the sum of emitted laser shots is displayed. It is also possible to display the applied energy. by pressing the READY LED. The cumulated energy in Joule is displayed for a short time (2-3s) and the device changes to STANDBY. The total energy is accumulated until the RESET button is pressed. It resets both the shot counter and the cumulated applied energy.

Note: Your service technician can activate or deactivate this feature.



Counter

RESET


6.2.7 DuoSpot Aiming Beam

The Q-LAS is equipped with a DuoSpot aiming beam as standard. The aiming beam is only visible in READY mode or while setting its brightness in STANDBY mode.


The two beams intersect in the focal plane to form a single spot. The laser exerts its effect on the focal plane. An exception is the focus shift (see chapter 6.3).

6.2.8 Head-Up Display

The Q-LAS has an integrated head-up display. As depicted below, the head-up display is able to display the selected power or the number of shots fired next to the aiming beam. This additional display allows you to keep track of the parameters during the procedure and avoid refocusing the treatment field. If you adjust the energy of the device, the head-up display shows the selected parameter in the upper right field of view for a few seconds. If a laser shot is fired, the number of emitted shots is displayed in the upper right field of view. You can adjust the brightness of the displayed values, using the rotating wheel at the upper right edge of the head-up display.

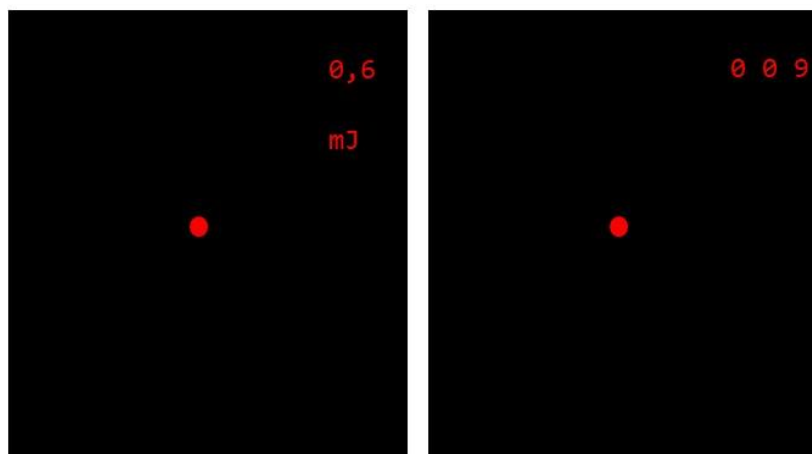

6.3 Adjustments at the laser head

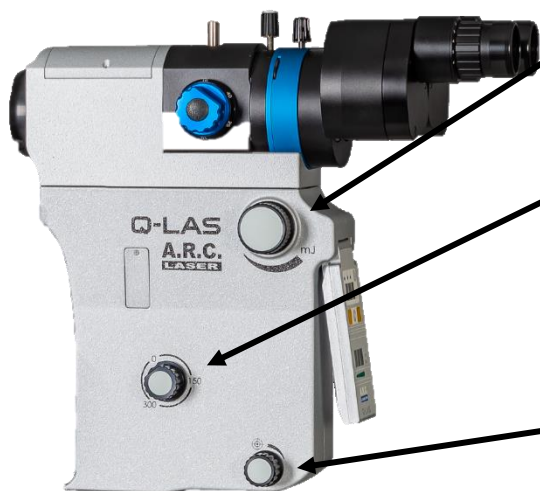

Power Selection

Use the energy knob to increase or decrease the energy emitted during one shot of the laser.

Focus Shift

The laser focus will be moved posterior:

0: Area of laser effect is ~30 µm behind the focal plane
150: Area of laser effect is 150 µm behind the focal plane
300: Area of laser effect is 300 µm behind the focal plane

Aiming Beam

Turn the knob to increase or decrease the brightness of the aiming beam.

6.4 Trigger Button at the Joystick

6.4.1 Standard version with rotating prism

To trigger the laser, press the trigger on the elevator (joystick) or optionally the foot switch. To do this, the laser must be switched to READY mode and the lighting arm must be rotated approx. 12° to the right or left. An acoustic signal and the display (POS) warn if the arm is not fully rotated out.



The laser shot can only be triggered when the arm is fully rotated out.

When a laser shot is fired, a warning signal sounds simultaneously and the LED next to the joystick lights up.

6.4.2 Version with lowering prism

To trigger the laser, press the trigger on the elevation (joystick) or optionally the foot switch. The laser must be switched to the READY mode for this. The prism must then be lowered by turning the adjustment wheel on the prism. An acoustic signal and the display indicator (POS) will be noticeable to warn the operator if the prism is not lowered.

Alternatively, the procedure described in Chapter 6.4.1 can be used.

When a laser shot is fired, a warning signal sounds simultaneously and the LED next to the joystick lights up.



6.5 Treatment

Position the patient on the chin rest. Make sure that the patient only touches the application parts. The application parts are the chin rest and the handles. Make sure that you and the patient do not touch any other parts at the same time.

Adjust the position of the slit lamp, illumination and focus so that the visible aiming beam points merge into one point.

Choose your initial values for energy and burst mode. It is recommended to start with a low energy level and then increase the energy during treatment.

Fix the patient's eye with a contact glass and focus the target beam into the eye to be treated.

Select the necessary beam intensity. The target point should be in focus with the microscope.

Press the READY button: the yellow LED flashes for 2 seconds and then changes to steady light. The treatment can start.



After setting the first treatment values, it is advisable to first set a laser pulse in the periphery of the retina, or the treatment area, in order to observe the laser effect. Under no circumstances should the first laser pulse be placed centrally or close to the macula.

ATTENTION

Any serious incident that occurs with this laser must be reported to the A.R.C. Laser GmbH and the responsible state authority.

6.6 Switch Off

- 1) Turn the key switch to position „OFF“(see chapter 4.10) to turn the laser off
- 2) Remove the key.
- 3) Switch off the unit with the main switch at the table.

ATTENTION

Do not forget to turn off the system. It should be covered with a dust cover to avoid dust deposits on the optics.

6.7 Laser-Stop

Restart the device after pressing the laser stop (red button):

- 1) Turn the key switch to position „OFF“
- 2) Pull out the red knob on the right side
- 3) Start the device again

The counter is reset automatically. The energy display shows the selected energy value.

ATTENTION

The **Laser-Stop** should only be used in emergencies.
It is located next to the key switch.

7 Specifications

7.1 General

Model Q-LAS	Q-switch Nd:YAG Laser
Cooling	passive with ambient air
Weight	44 kg with table and slit lamp
Dimensions	H 130 cm/ B 86 cm/ T 57 cm with slit lamp, height adjustable table 72-94 cm
Control Option	7-segment display

7.2 Laser Data

Wavelength	Energy	Tolerance
1064 nm	30 mJ max. (BurstMode 3) Settings (BurstMode 1): 0.1 mJ Steps from 0.5 to 10 mJ (up to 20 mJ – BurstMode 2) BurstMode: 1,2 or 3 Pulses	±20% of the shown energy according to IEC 60601-2-22 extended tolerance range 20 to 30 % of the shown energy according to IEC 60601-2-22 Abs.201.13.2.101

Output	Single pulse, repetition rate 3 Hz, 2.4 Hz or 1.8 Hz
Pulse duration	< 5ns
Aiming beam	635 nm red < 1mW, variable brightness
Power transmission	coupled with a slit lamp

7.3 Electrical Connection Data

Power supply connection values	100-240 V, 50Hz/60 Hz, 250 VA
--------------------------------	-------------------------------

7.4 Classification

Laser class Laser Beam (Classification EN 60825-1)	3B
Laser class Aiming Beam (Classification EN 60825-1)	2
Classification according to MDR	IIb Rule 9
Electrical protection class (Classification IEC 60601-1)	II
Certification	CE 0123

8 Service

8.1 Introduction

The device was designed, developed and tested according to the latest technical knowledge. We have set the product life to 7 years. In addition, the availability of spare parts is guaranteed by us within a period of 10 years. However, in order to ensure that everything works properly, we have made it possible for you to carry out a visual check of the status indicators from the outside.

ATTENTION

There is no need for the laser user to perform routine or service work within the laser system. All adjustments and calibrations that require the protective housing must be carried out by trained service personnel. This also includes cleaning and cleaning the optics within the laser.

8.2 Safety Check (STK)

Once within 24 months, the laser should be subjected to a safety check (STK) by trained personnel.

Please note that a medical device book and a regular STK safety check are not required in every country. Note local requirements and laws.

Scope of Safety Check (STK)

Visual Check

- Laser marking (laser class, max. Power, wavelength)
- Information signs/warning signs; properly and completely attached
- Instructions for use/medical device book
- Condition of the supply line
- Condition of the goggles/protective device
- Overall condition

Functionality Check

- Footswitch/joystick
- Beam guidance system/coupling/decoupling/pilot laser
- Check operating foil (touch screen)

Check the necessary monitoring/safety display and signaling unit

- Laser safety glasses
- Emission controls (acoustic, visual)
- Power meter (comparison internal, external)
- Key switch
- Laser Stop switch (check for function)
- Interlock device (check for function)

Electrical Safety Check

- According to IEC 62353 or DIN VDE 0751, Part 1

Output power measurements

- Check the set power with an external power meter (permissible tolerance + 20%)

ATTENTION

If one or more safety-related points are objected to after the safety-related inspection (STK), the device may no longer be operated.

8.3 Care and Maintenance by the User

The following care instructions can be carried out by the user. These serve to make your work easier. For cleaning, the system must be disconnected from the mains. Always use a damp, but never wet, soft cloth for cleaning and disinfecting.

For cleaning and disinfection, the power supply must be disconnected. First clean the device with clear water and neutral detergent to remove coarse and visible contamination. Make sure that no moisture gets into the device.

Furthermore, wipe disinfection is possible. When choosing the disinfectant, pay attention to the following:

- according to the manufacturer of the disinfectant, the disinfectant should be suitable for non-invasive medical devices
- according to the manufacturer of the disinfectant, the disinfectant should be suitable for wipe disinfection of surfaces
- the disinfectant should be based on alcohol and/or quaternary compounds
- the disinfectant should be suitable for lacquers
- aldehyde-free disinfectant (recommended)

When doing wiping disinfection, the manufacturer's instructions for the disinfectant must be observed. Following possible disinfectant that meet the above mentioned requirements are listed:

Manufacturer	Possible disinfectants*
BODE Chemie GmbH	Mikrobac forte, Bacillol 30 Foam
ANTISEPTICA	Acrylan, Biguacid Liquid
Schülke & Mayr GmbH	acryl-des, antifact AF (N)
Dr. Schumacher GmbH	CLEANISEPT, Descosept AF
Ecolab	Incidin Foam, Incidin Pro
Dr. Weigert	neoform MED AF, neoform MED rapid

*This list is not exhaustive.

ATTENTION

When cleaning, the device must always be switched off and disconnected from the mains. Wet wipes should be avoided in any case. Exposure to water can lead to defects.

8.3.1 Inspection and cleaning of external optics

Check the accessible surfaces of the optics for possible contamination. The surfaces of the slit lamp can be cleaned with a soft cloth and distilled water or a mixture of distilled water and approx. 10% alcohol.

8.3.2 Inspection and cleaning of internal Slit lamp optics

Regularly check all optics on the slit lamp - especially the plug-in mirror and the lens above it - for contamination.

A contaminated prism (in PCL Z slit lamp) could lead to a reduction of energy when the laser beam passes through.

Remedial Action:

- 1) Use non-alcoholic, not soaked cotton swabs and try to clean the outside of the lens with distilled water. If this does not lead to success, please continue with step number 2.
- 2) Take a few drops of methanol and gently drizzle the end of the cotton swab, then do as in step number 1.

ATTENTION

Always wipe the optical surface in one direction. Never go back in the other direction with the same cotton swab. Any particles that have already settled in the cotton swab would scratch the optical surface again when moving back and thus cause irreparable damage to the optics.

8.4 Slit lamp

You can find detailed information on the care and maintenance of the slit lamp in the accompanying manual. It is essential to follow these instructions.

8.4.1 Dust Cover

After each use, the dust cover and the red display protective cover should be replaced to keep the surfaces clean.

8.4.2 Test the Slit lamp focus

- 1) Use the focus rod and insert it into the hole in the slit lamp axis. Follow the slit lamp manual that comes with the device.
- 2) Turn on the laser. Make sure that the laser is not triggered during this time. The laser must be in the READY mode. Center the aiming beam in the middle of the visual field on the focus rod.
- 3) Check that the aiming beam is circular and symmetrical.
- 4) A disturbed target beam can be an indication of a faulty system. Checking the aiming beam for indications of the integrity of the system.



ATTENTION

Press the eyepieces gently into the holder so that they are held completely in the binocular piece. Adjust each eyepiece to your satisfaction so that you can see a sharp imager of the focus rod through each eyepiece.

8.5 Troubleshooting

Possible Error		Messages and solutions	
F00	Watchdog error	Device is hanging during system check	Hardware Error Please restart the system. If the error occurs regular, contact your local sales & service-representative of A.R.C. Laser
F01	Reduced output -20%	Permanently checked by test shots during the system start-up and during the operation	Hardware error Please restart the system
F02	Higher output +30%	Permanently checked by test shots during the system start-up and during the operation	Hardware error Please restart the system
F04	Optokoppler Safety Shutter Time Out	Permanently checked by test shots during the system start-up and during the ready phase	Please contact your local sales & service-representative of A.R.C. Laser
F06 F07 F08	Foot switch short circuit – pre / post	Checked during the system start	Please contact your local sales & service-representative of A.R.C. Laser
F11	Unable to set power	Checked during the system start	Please contact your local sales & service-representative of A.R.C. Laser
F21	Error READY button	Checked during the system start	Please contact your local sales & service-representative of A.R.C. Laser
F38	DA converter	Checked during the system start	Please contact your local sales & service-representative of A.R.C. Laser
F44	Error pulse monitoring and internal temperature	Checked during the system start	Please contact your local sales & service-representative of A.R.C. Laser
F47	Performance check, Checksum error	Checked during the system start	Please contact your local sales & service-representative of A.R.C. Laser
F50	I2C-error	Checked during the system start	Please contact your local sales & service-representative of A.R.C. Laser
F51	Keys error	User interface error at the front panel	Please contact your local sales & service-representative of A.R.C. Laser
F58	Diode current too high	Checking the maximum diode current	Please contact your local sales & service-representative of A.R.C. Laser
F62	Energy test error or burst mode voltage monitoring error	Checking if the energy or the burst voltage monitoring is function properly (d)	Please contact your local sales & service-representative of A.R.C. Laser
F63	Temperature/ humidity error	Checking the temperature, humidity respectively	Decrease the temperature, humidity respectively.
F65	Connection-loss of the Head-up display	Checking the connection/signal of the Head-up display during the energy settings	Please contact your local sales & service-representative of A.R.C. Laser
FOC	Position of the focus shift dial	Position of the focus shift is constantly checked	Check the position of the focus shift. It should be locked.
IL1	Internal Temperature	The internal temperature is out of permitted range. The laser should be allowed to cool down.	Please press the Reset-Button to delete the message.
IL2	Time Delay Trigger	Pressing the trigger only half way down or too slowly may cause this error message.	Please press the Reset-Button to delete the message. You immediately can proceed working afterwards.
IL3	Internal power monitoring	Error when changing the energy while the device changes from STANDBY to READY	Device goes to STANDBY. Press the counter reset button to clear the display.

Possible Error		Messages and solutions	
POS	Position of the arm	The position of the lighting arm is constantly checked	Swing out the arm

8.6 Disposal

The relevant, locally applicable laws and regulations must be observed during disposal. Under no circumstances should the device be disposed of with domestic waste.

A.R.C. Laser GmbH is happy to assist with the disposal of the system; Costs for the proper return of the laser to A.R.C. Laser GmbH is the responsibility of the customer. Please contact our service department.



9 Customer Service

9.1 Warranty Information

A.R.C. Laser GmbH grants you a two-year guarantee. Parts that have a defect will be replaced free of charge within two years. All add-on and purchased parts are exempt from this guarantee. Our guarantee extends to the repair or replacement of defective parts. However, we reserve the right to renew entire assemblies and adapt them to technical progress.

Repairs by third parties or changes to the device will void the warranty. The use of other parts that have not been accepted with the device or obtained from other suppliers will also void the warranty. The attachment of parts or assemblies or other changes to the device requires the express written confirmation by A.R.C. Laser GmbH.

9.2 Warranty, Shipment, Packing

A warranty claim for defective parts, malfunction or damage to the housing of the device must be submitted to A.R.C. Laser GmbH within 24 hours. Parts that are returned during the warranty period (at the express request of A.R.C. Laser GmbH) must be confirmed in writing by A.R.C. Laser GmbH. Detailed packaging instructions and information on how to return the device are provided by A.R.C. Laser GmbH. The return must be insured and paid for by shipper. The insurance and transportation costs are not covered by A.R.C. Laser GmbH. The choice of the return is made by the A.R.C. Laser GmbH communicated to the customer. Changes and amendments in the carrier or the shipping method can lead to delays in transport and processing. All components to be changed under the warranty claim are manufactured by A.R.C. Laser GmbH renewed free of charge within the guarantee period. We reserve the right to make changes to the design of the device - if it appears necessary - to increase the safety or the functionality of the device. The responsibility for the design as well as for changes in the device lies solely with A.R.C. Laser GmbH. Changes will be communicated to the customer and accordingly carried out at A.R.C. Laser GmbH.

9.3 Sales and Service Information

For sales and service information, please contact A.R.C. Laser GmbH or our local distributor.

10 Guidelines and Manufacturers Declaration

The Q-LAS is intended for use in the professional environment of a healthcare facility. The Q-LAS does not have essential performance characteristics.

10.1 Electromagnetic Emissions

The laser is intended for use in an environment as specified below. The customer or user of the laser should ensure that it is operated in such an environment.		
Immunity tests	Compliance	Electromagnetic environment - guideline
RF-Emissions CISPR 11 150 kHz – 30 MHz	EN 55011 Group 1/Class B	The laser uses RF energy exclusively for its internal function. Hence, RF emission is very low and not likely to cause any interference in nearby electronic equipment
RF- Emissions CISPR 11 30 MHz - GHz	EN 55011 Group 1/Class B	The laser is only suitable for the environment in professional healthcare facilities
Harmonic emission	IEC 61000-3-2 Class A	
Voltage fluctuations/flicker	IEC 61000-3-3	


10.2 Electromagnetic Immunity (1)

The laser is intended for use in the electromagnetic environment specified below. The customer or the user of the laser should ensure that it is used in such an environment.			
Immunity tests	IEC 60601-Test level	Compliance level	Electromagnetic environment – guidelines
Electrostatic discharge (ESD)	± 2 kV, ± 4 kV, ± 6, ± 8 contact discharge; ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air discharge	IEC 61000-4-2	Floors should be made of wood, ceramic or stone. If the floor is covered with a synthetic material, the relative air humidity should be at least 30%.
Radiated RF electromagnetic fields	3 V/m 80 MHz to 2.7 GHz 80 % AM at 1 kHz	IEC 61000-4-3	Only the voltage of a typical professional healthcare facility may be used..
Proximity fields from RF wireless communications equipment	3 V/m 80 MHz to 2,7 GHz (see Table 10.4)	IEC 61000-4-3	Only the voltage of a typical professional healthcare facility may be used..
Rated power frequency magnetic fields	30 A/m 50 Hz or 60 Hz	IEC 61000-4-8	Magnetic fields at the grid frequency that are usually available in domestic or clinical areas may be used.
Electrical fast transients /burst	± 2 kV for power lines ± 1 kV for IO-lines 100 kHz repetition frequency	IEC 61000-4-4	Only the voltage of a typical professional healthcare facility may be used..
Surge voltages (Surges), Line against line	± 0.5 kV, ± 1 kV	IEC 61000-4-5	Only the voltage of a typical professional healthcare facility may be used.
Surge voltages (Surges), Line against grounding	± 0.5 kV, ± 1 kV, ± 2 kV	IEC 61000-4-5	Only the voltage of a typical professional healthcare facility may be used.
Conducted disturbances induced by RF fields	3 V 0.15 MHz to 80 MHz 6 V in ISM-frequency bands between 0.15 MHz and 80 MHz 80% AM at 1 kHz	IEC 61000-4-6	Only the voltage of a typical professional healthcare facility may be used.
Voltage dips	0 % UT; ½ cycle at 0, 45, 90, 135, 180, 225, 270 and 315 level 0 % UT; 1 cycle at 0 and 180 level And 70 % UT; 25/30 cycles at 0 and 180 level	IEC 61000-4-11	Floors should be made of wood, ceramic or stone. If the floor is covered with a synthetic material, the relative air humidity should be at least 30%.

Voltage interruption	0% UT; 250/300 cycles at 0 and 180 level	IEC 61000-4-11	<p>Quality of the supply voltage that is usually available in domestic or clinical areas must be used.</p> <p>It is recommended to use an uninterruptible power supply.</p>
NOTE: UT is the AC mains voltage prior to application of the test level.			

10.3 Electromagnetic Immunity (2)

The laser is intended for use in the electromagnetic environment specified below. The customer or the user of the laser should ensure that it is used in such an environment.

Immunity tests	IEC 60601-Test level	EMV standard	Electromagnetic environment – guidelines
Conducted disturbance variables, induced by high-frequency fields	3 V/m 150 kHz to 80 MHz	IEC 61000-4-6	<p>Portable and mobile RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the CLASSIC 514, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.</p> <p>The field strength of stationary radio transmitters is, as determined by an electromagnetic site survey, at all frequencies smaller than the compliance level.</p>
Radiofrequency electromagnetic fields in the immediate vicinity of wireless communication devices	3 V/m 80 MHz to 2.7 GHz	IEC 61000-4-3	<p>Interference may occur in the vicinity of equipment marked with the following symbol:</p>  <p>If higher IMMUNITY TEST LEVELS than those specified in Table 9 (IEC 60601-1-2) are used, the minimum separation distance may be lowered. Lower minimum separation distances shall be calculated using the equation specified in 8.10 (IEC 60601-1-2).</p>

10.4 Recommended separation distances between portable and mobile RF telecommunications equipment and the laser

The laser is intended for use in an electromagnetic environment in which the RF disturbances are controlled. The customer or the user of the laser can help to avoid electromagnetic interference by maintaining the minimum distance between portable and mobile HF telecommunication devices (transmitters) and the laser - depending on the output power of the communication device, as stated below.

Test frequency	Frequency band	Radio service	Modulation	Maximum performance	Distance	Immunity test level
MHz	MHz			W	m	V/m
385	380 to 390	TETRA 400	Pulse modulation 18Hz	1.8	0.3	27
450	430 to 470	GMRS 460 FRS 460	FM ± 5 kHz Hub 1 kHz Sinus	2	0.3	28
710 745 780	704 to 787	LTE Band 13,17	Pulse modulation 217 Hz	0.2	0.3	9
810 870 930	800 to 960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation 18 Hz	2	0.3	28
1720 1845 1970	1700 to 1990	GSM 1800, CDMA 1900, GSM 1900, DECT, LTE Band 1,3,4,25, UMTS	Pulse modulation 217 Hz	2	0.3	28
2450	2400 to 2570	Bluetooth, WLAN 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation 217 Hz	2	0.3	28
5240 5500 5785	5100 to 5800	WLAN 801.11 a/n	Pulse modulation 217 Hz	0.2	0.3	9

Notes:

:



A.R.C. LASER GmbH

Bessemerstr. 14
D-90411 Nürnberg

Tel: +49 911 21779-0
Fax: +49 911 21779-99

www.arclaser.de